MICROACTIVE TRISTAR® II PLUS

SURFACE AREA AND POROSITY ANALYZER



micromeritics®

CALCULATIONS

Oct 2020

CONTACT US

Micromeritics Instrument Corporation

4356 Communications Drive Norcross, GA / USA / 30093-2901 Phone: 1-770-662-3636 Fax: 1-770-662-3696 www.Micromeritics.com

Instrument Service or Repair

Phone: 1-770-662-3666 International — contact your local distributor or call 1-770-662-3666 Service.Helpdesk@Micromeritics.com

Micromeritics Learning Center

Phone: 1-770-662-3607 www.Micro.edu

CORPORATE PROFILE

Micromeritics Instrument Corporation is a leading global provider of solutions for material characterization with bestin-class instrumentation and application expertise in five core areas: density; surface area and porosity; particle size and shape; powder characterization; and catalyst characterization and process development.

Founded in 1962, the company is headquartered in Norcross, Georgia, USA and has more than 400 employees worldwide. With a fully integrated operation that extends from a world class scientific knowledge base through to inhouse manufacture, Micromeritics delivers an extensive range of high-performance products for academic research and industrial problem-solving.

Under its premium brand Particulate Systems, Micromeritics discovers and commercializes unique and innovative material characterization technologies that are complementary to core product lines.

The company's holistic, customer-centric approach also extends to a cost-efficient contract testing laboratory – the Particle Testing Authority (PTA). The strategic acquisitions of Freeman Technology Ltd and Process Integral Development S.L. (PID Eng & Tech) reflect an ongoing commitment to optimized, integrated solutions in the industrially vital areas of powders and catalysis.

Freeman Technology (Tewkesbury, UK) brings market-leading powder characterization technology to Micromeritics' existing portfolio of particle characterization techniques. The result is a suite of products that directly supports efforts to understand and engineer particle properties to meet powder performance targets. With over 15 years of experience in powder testing, Freeman Technology specializes in systems for measuring the flow properties of powders. In combination with detailed application know-how these systems deliver unrivalled insight into powder behavior supporting development, formulation, scale-up, processing and manufacture across a wide range of industrial sectors.

PID Eng & Tech (Madrid, Spain) complements Micromeritics' renowned offering for catalyst characterization with technology for the measurement and optimization of catalytic activity, with a product range that extends to both standard and bespoke pilot scale equipment. Launched in 2003, PID Eng & Tech is a leading provider of automated, modular microreactor systems for the detailed investigation of reaction kinetics and yield. These products are supported by a highly skilled multidisciplinary team of engineers with in-depth expertise in the design, construction and operation of laboratory units and process scale-up.

The Particle Testing Authority (PTA) provides material characterization services for fine powders and solid materials using Micromeritics' instrumentation alongside complementary solutions from other vendors. With the certification and expertise to operate across a wide range of industries the PTA offering runs from single sample analysis to complex method development, method validation, new product assessment, and the analytical support required for large scale manufacturing projects. An experienced, highly trained team of scientists, engineers, and lab technicians works closely with every client to ensure that all analytical requirements are rapidly and responsively addressed.

Micromeritics has a strong global network with offices across the Americas, Asia, and Europe complemented by a dedicated team of distributors in additional locations. This ensures that local, knowledgeable support is available for every customer, in academia or industry. Micromeritics works across a truly diverse range of industries from oil processing, petrochemicals and catalysts, to food and pharmaceuticals, and at the forefront of characterization technology for next generation materials such as graphene, metal-organic-frameworks, nanocatalysts, and zeolites. Engineering solutions that work optimally for every user is a defining characteristic of the company.

mi micromeritics[®]

Table of Contents

Contact Us	<i>ii</i>
Corporate Profile	
Alpha-S Method	1
BET Surface Area	2
BJH Pore Volume and Area Distribution	3
DFT (Density Functional Theory)	12
Dollimore-Heal Pore Volume and Area Distribution	15
Dubinin-Astakhov	17
Dubinin-Radushkevich	21
Equilibration	23
f-Ratio Method	24
Freundlich Isotherm	25
Heat of Adsorption for Physical Adsorption	26
Horvath-Kawazoe	27
Langmuir Surface Area	37
MP-Method	39
Quantity Adsorbed	41

mi micromeritics[®]

Relative Pressure	43
Saturation Pressure	44
SPC Report Variables	45
Summary Report	47
t-Plot	49
Temkin Isotherm	50
Thermal Transpiration Correction	51
Thickness Curve	52

Blank Page

ALPHA-S METHOD

The alpha-S curve is calculated from the reference isotherm by dividing each quantity adsorbed by the quantity adsorbed at 0.4 relative pressure.

$$lpha_i = rac{Q_i}{Q_{0.4}}$$

where $Q_{0.4}$ is found by linear interpolation.

A least-squares analysis fit is performed on the $(\alpha_i, Q_{ads,i})$ data pairs where α_i is the independent variable and $Q_{ads,i}$ is the dependent variable. The following are calculated:

- a. Slope (S cm³/g STP)
- b. Y-intercept (Q₀ cm³/g STP)
- c. Error of the slope (cm^3/g STP)
- d. Error of the Y-intercept (cm^3/g STP)
- e. Correlation coefficient

Surface area is calculated as:

$$A_s = rac{A_{
m ref}S}{Q_{0.4}}$$

Pore size is calculated as:

 $\frac{Q_0V_{\rm mol}}{22414{\rm cm}^3{\rm STP}}$

where V_{mol} is liquid molar volume from the fluid property information.

BET SURFACE AREA

The BET¹) transformation is calculated as:

$$y=rac{1}{Q(
ho^0/
ho-1)}$$

A least-squares fit is performed on the $P_{\rm rel}$, y. The following are calculated:

- a. Slope (S cm³/g STP)
- b. Y-intercept (y₀ cm³/g STP))
- c. Uncertainty of the slope $(u(s) \text{ cm}^3/\text{g STP}))$
- d. Uncertainty of the Y-intercept $(u(y_0) \text{ cm}^3/\text{g STP}))$
- e. Correlation coefficient

Using the results of the above calculations, the following can be calculated:

BET Surface Area (m²/g):

$$A_s = rac{A_m N_A}{V_m (s+y_0)} imes 10^{-18} m^2 / nm^2$$

BET C value:

$$C = rac{s}{y_0} + 1$$

Quantity in the Monolayer (cm³/g STP):

$$Q_m = rac{1}{Cy_0} = rac{1}{s+y_0}$$

Error of the BET Surface Area (m^2/g) :

$$u(A_s) = rac{\sqrt{u^2(s) + u^2(y_0)}}{s + y_0}$$

¹) Brunauer, S.; Emmett, P.H.; and Teller, E., *J. Am. Chem. Soc.* 60, 309 (1938).

mi micromeritics[®]

BJH PORE VOLUME AND AREA DISTRIBUTION

For adsorption data, the relative pressure and quantity adsorbed data point pairs collected during an analysis must be arranged in reverse order from which the points were collected during analysis. All calculations are performed based on a desorption model, regardless of whether adsorption or desorption data are being used.

The data used in these calculations must be in order of strictly decreasing numerical value. Points which do not meet this criterion are omitted. The remaining data set is composed of relative pressure (*P*), quantity adsorbed (*Q*) pairs from (P_1 , Q_1) to (P, Q_n) where ($P_n = 0$, $Q_n = 0$) is assumed as a final point. Each data pair represents an interval boundary (or desorption step boundary) for intervals i=1 to i=n-1 where n = total number of (*P*, *Q*) pairs.

Generally, the desorption branch of an isotherm is used to relate the amount of adsorbate lost in a desorption step to the average size of pores emptied in the step. A pore loses its condensed liquid adsorbate, known as the core of the pore, at a particular relative pressure related to the core radius by the Kelvin¹) equation. After the core has evaporated, a layer of adsorbate remains on the wall of the pore. The thickness of this layer is calculated for a particular relative pressure from the thickness equation. This layer becomes thinner with successive decreases in pressure, so that the measured quantity of gas desorbed in a step is composed of a quantity equivalent to the liquid cores evaporated in that step plus the quantity desorbed from the pore walls of pores whose cores have been evaporated in that and previous steps. Barrett, Joyner, and Halenda²) developed the method (known as the BJH method) which incorporates these ideas. The algorithm used is an implementation of the BJH method.


EXPLANATION OF TERMS

A pore filled with condensed liquid nitrogen has three zones:

- **Core.** The *core* evaporates all at once when the critical pressure for that radius is reached; the relationship between the core radius and the critical pressure is defined by the Kelvin equation.
- Adsorbed layer. The adsorbed layer composed of adsorbed gas that is stripped off a bit at a
 time with each pressure step; the relationship between the thickness of the layer and the relative
 pressure is defined by the thickness equation.
- Walls of the cylindrical pore. The *walls of the cylindrical pore* the diameter of the empty pore is required to determine the pore volume and pore area. End area is neglected.

¹) Kelvin, J. (published under the name of Sir William Thomson), *Phil. Mag.* 42, 448-452 (1871).

²⁾ Barrett, E.P.; Joyner, L.S.; and Halenda, P.P., J. Am. Chem. Soc. 73, 373-380 (1951).

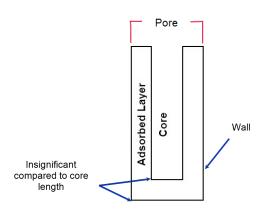

CALCULATIONS

The quantities adsorbed (Qa) are converted to the liquid equivalent volumes (VI, cm³/g)

$$Vl_i = rac{Q_i V_{
m mol}}{22414 {
m cm}^3 {
m STP}}$$

where V_{mol} is the liquid molar volume from the fluid property information.

The relative pressure (P_i) is assumed to be close to unity so that substantially all the pores in the sample are filled.

The corresponding Kelvin core radius is calculated. Only pores smaller than this size will be included:

$$Rc_i = rac{-A}{(l+F)\ln(P_i)}$$

where

- A = adsorbate property factor (from the *BJH Adsorptive Options* window)
- F = fraction of pores open at both ends (from the BJH Adsorption Report Options window or the BJH Desorption Report Options window); assumed to be zero for desorption

Rc = Kelvin radius (Å) of core

This radius will be adjusted for the thickness of the adsorbed layer during subsequent calculation steps.

mi micromeritics[®]

The following calculations (a-c) are made for each relative pressure interval based on the increment of volume desorbed during that interval:

where

- i = interval number, that is i=1 for the first interval from P1 to P2, and so on
- j = each previous interval during which new pores were found
- k = the total number of intervals in which new pores have been found. It is also the number of lines reported on the BJH table for collected data
- a. The thickness of the adsorbed layer at the end of the interval is calculated using the equation located in *Thickness Curve on page 52*.

For the last pressure interval from the lowest Pr_i to zero relative pressure, reference the calculations from the equations in <u>Thickness Curve on page 52</u>.

For the first pressure interval, there are no previously opened pores so the volume of liquid desorbed from walls of previously opened pores is zero ($Vd_1 = 0$), and the remainder of Step (a) is skipped.

The change in thickness of the wall layer due to desorption from previously opened pores is calculated as:

 $\Delta Tw = Tw_1 - Tw_{i+1}$

The annular cross-sectional area of the wall layer desorbed is calculated for all previously opened pores:

$$CSA_{
m j}=\piiggl[iggl(Rc_{
m j}=\Delta Twiggr)^2-Rc_{
m j}^2iggr]iggl(10^{-16}rac{{
m cm}^2}{{
m \AA}^2}iggr)$$

The total volume of gas desorbed from walls of previously opened pores is calculated:

 $Vd_{
m i} = rac{\sum}{
m j} (LP_j)(CSAa_{
m j})$ for all previously opened pores

where LP_i = length of previously opened pores as calculated in Step b(2).

- b. The physical processes occurring for this pressure interval are determined as:
 - 1. If Vd_i is greater than the current increment of volume desorbed $(Vl_i Vl_{i+1})_i$, desorption from walls only is occurring. Total surface of walls exposed thus far (cm²/g) is calculated as:

$$SA_{\mathbf{w}} = rac{\sum}{\mathrm{j}} \pi (LP_{\mathrm{j}}) (D_{\mathrm{avg,j}}) (rac{\mathrm{10^{-8} cm}}{\mathrm{\dot{A}}})$$

for all previously opened pores

where

 $D_{\text{avg,j}}$ = weighted average pore diameter calculated in Step b.2.

A new layer thickness (ΔTw) that will not overcompensate for the actual volume desorbed in this interval is calculated:

$$\Delta Tw = rac{(Vl_{ ext{i}}-Vl_{ ext{i}+ ext{i}})(10^8rac{\dot{ ext{A}}}{ ext{cm}})}{SAw_{ ext{i}}}$$

Since no cores are evaporated in this pressure interval, no new pores are revealed. Thus no ending Kelvin radius and average pore diameter are calculated for this interval. Note that this means the report may have fewer tabulated intervals on the collected data report than experimental pressure intervals.

2. If Vd_i is less than the volume increment desorbed during this interval $(Vl_i - Vl_{i+1})$, the remaining volume is due to new pores with core evaporation taking place in this interval. *K*, the number of intervals with new pores exposed, is increased by 1. (For the interval from the lowest Pr_1 to zero relative pressure, no new pore volume is calculated and the rest of Step b is skipped.)

The volume desorbed from newly opened pores in this interval is calculated as:

$$Vc_{\mathrm{i}} = (Vl_{\mathrm{i}} - Vl_{\mathrm{i+1}}) - Vd_{\mathrm{i}}$$

The Kelvin radius for the end of the interval is calculated as:

$$Rc_{k+1} = rac{-A}{(1+F)\ln(P_{i+1})}$$

All new pores opened in this interval are represented by one pore having a lengthweighted average pore diameter and a corresponding length sufficient to account for the required volume of adsorbate. The weighted average pore diameter is calculated as:

$$D_{\text{avg},\text{k}} = rac{2(Rc_{\text{k}}+Rc_{\text{k}+1})(Rc_{\text{k}})(Rc_{\text{k}+1})}{Rc_{\text{k}}^{2}+Rc_{\text{k}+1}^{2}}$$

mi micromeritics[®]

 $D_{avg,k}$ is the diameter of a pore which would have a surface area that is the average of the areas for pores radius Rc_k and Rc_{k+1} , if its length was the mean of the lengths at those radii.

The relative pressure corresponding to $D_{avg,k}$ is calculated as:

$$P_{\mathrm{avg,k}} = \exp\left[rac{-2A}{(1+F)(D_{\mathrm{avg,k}})}
ight]$$

The thickness of the adsorbed layer at this pressure is calculated as:

$$Tw_{\mathrm{avg,k}} = HP1 iggl[rac{HP2}{ln(P_{\mathrm{avg,k}})}iggr]^{HP3}$$

The decrease in thickness of the wall layer by desorption from the walls of new pores during the lower portion of the pressure interval is calculated as:

$$\Delta T d = T w_{ ext{avg,k}} - T w_{ ext{i+1}}$$

The cross-sectional area of the newly opened pores is calculated as:

$$CSAc_{ extsf{k}} = \left[rac{D_{ extsf{avg,k}}}{2} + \Delta Td
ight]^2 igg(rac{10^{-16} extsf{cm}^2}{ extsf{A}^2}igg)$$

The length of the newly opened pores is calculated as:

$$LP_{\mathbf{k}} = rac{Vc_{\mathbf{i}}}{CSAc_{\mathbf{k}}}$$

Pore diameters and radii are adjusted for the change in thickness of the adsorbed wall layer during this interval. If new pores were opened during this interval, the average diameter is adjusted by the change in layer thickness during the second portion of the desorption interval as:

$$D_{\mathrm{avg,k,new}} = D_{\mathrm{avg,k,old}} + 2(\Delta T d)$$

mi micromeritics[®]

The layer thickness change during the whole interval is added to diameters of previously opened pores as:

$$D_{\mathrm{avg,k,new}} = D_{\mathrm{avg,k,old}} + 2(\Delta T dw)$$

(not including Davg,k)

The layer thickness change desorbed during this interval also is added to the radii corresponding to the ends of the pressure intervals as:

$$Rc_{ ext{j,new}} = Rc_{ ext{j,old}} + \Delta Tw$$

for all except Rck+1.

Steps *a* to *c* are repeated for each pressure interval.

After the above calculations have been performed, the diameters corresponding to the ends of the intervals are calculated as:

$$Dp_{
m j}=2(Rc_{
m j})$$

for all Rc_i . including Rc_{k+1} .

The remaining calculations are based on Dp_i , $D_{avg,i}$, and LP_i . These calculations are only done for $D_{avg,i}$ values that fall between the Minimum BJH diameter and the Maximum BJH diameter specified by the operator on the *BJH Adsorption Report Options* window or the *BJH Desorption Report Options* window.

(1) Incremental Pore Volume (Vp_i, cm³/g):

$$Vp_{\mathrm{i}}=\pi(Lp_{\mathrm{i}})iggl[rac{D_{\mathrm{avg,i}}}{2}iggr]^2iggl[rac{10^{16}\mathrm{cm}^2}{\mathrm{\AA}^2}iggr]$$

(2) Cumulative Pore Volume ($VP_{cum,i}$, cm³/g):

$$VP_{ ext{cum,i}} = rac{\sum}{j} VP_{ ext{j}} ext{ for } (J \leq 1)$$

(3) Incremental Surface Area SA_{cum.i},m²g-A):

$$SA_{
m i}=\piigg(LP_{
m i}igg)igg(rac{10^{-2}{
m m}}{
m cm}igg)igg(D_{
m avg,i}igg)igg(rac{10^{-10}{
m m}}{
m \dot{A}}igg)$$

mi micromeritics®

(4) Cumulative Surface Area SA_{cum,i},m²g):

$$SA_{ ext{cum},10} = \sum SA_{ ext{j}} ext{for} J \leq 1$$

(5) dV/dD pore volume (dV/dD_i, $cm^3/g-A$):

$$rac{dV}{dD_{\mathrm{i}}} = rac{VP_{\mathrm{i}}}{Dp_{i}-Dp_{\mathrm{i+1}}}$$

(6) $dV/d \log (D)$ pore volume ($dV/d \log(D)_i$, cm³g):

$$rac{dDv}{d\log D_{
m i}} = rac{VP_{
m i}}{\log\left(rac{Dp_{
m i}}{Dp_{
m i+1}}
ight)}$$

(7) dA/dD pore area (dA/dD_i, m^2/g -A):

$$rac{dA}{dD\mathrm{i}} = rac{SA_\mathrm{i}}{Dp_\mathrm{i} - Dp_\mathrm{i+1}}$$

(8) *dA/d log(D)* pore area [*dA/d log(D)i*, m²/g]:

$$rac{dA}{d\log D_{\mathrm{i}}} = rac{SA_{i}}{\log\left(rac{D_{p_{\mathrm{i}}}}{D_{p_{\mathrm{i}+1}}}
ight)}$$

For fixed pore size tables (if selected), the following calculations are performed

(1) Average Fixed Pore Size (DF_{avg,j},A):

$$DF_{\mathrm{avg,j}} = rac{Dp_{F_{\mathrm{j}}} + D_{P_{F_{\mathrm{j+1}}}}}{2}$$

calculated for all intervals in the fixed pore size table.

For the intervals with between the Minimum BJH diameter and the Maximum BJH diameter.

(2) Cumulative Pore volume ($VpF_{cum,i}, cm^3/g$):

 $V_P F_{\mathrm{cum,i}} = \mathrm{INTERP}(D_P F_{\mathrm{i+1}})$

where INTERP(x) is the value interpolated from the function $X=dP_{j+1}$ and $Y=VP_{cum,i}$, using an AKIMA semi-spline interpolation.

mi micromeritics[®]

(3) Incremental Pore Volume (VpF_i, cum³/g):

$$V_P F_i = V_P F_{ ext{cum,i}} - V_P F_{ ext{cum,i-1}}$$

where VpF_{cum,0}=0

(4) Cumulative Surface Area (SAF_{cum,i}, m²/g):

$$SAF_{ ext{cum}, ext{i}} = ext{INTERP}(D_PF_{ ext{i}+1})$$

where INTERP(x) is the value interpolated from the function $X = Dp_{j+1}$ and $Y=SA_{cum,j}$.

(5) Incremental Surface Area (SAF_i, m²/g):

$$SAF_{
m i} = SAF_{
m cum,i} - SAF_{
m cum,i-1}$$

where $SAF_{cum.0} = 0$.

(6) dV/dD pore volume ($dV/dDpF_i$, cm³/g-A):

$$rac{dV}{dDpF_{
m i}} = {
m INTERP}(DpF_{
m i+1})$$

where INTERP(x) is the value interpolated from the function $X = D_{avg,j}$ and $Y=dV/dD_{j}$

(7) dV/d log(D) pore volume $[dV/d log(DpF_i), cm^3/g]$:

$$rac{dV}{d\log(DpF_{
m i})} = {
m INTERP}(DpF_{
m i+1})$$

where INTERP(x) is the value interpolated from the function $X = D_{avg,j}$ and Y=dV/d log(D)_i

(8) dA/dD pore area ($dA/dDpF_i$, m²/g-A):

$$rac{dA}{dDpF_{
m i}} = {
m INTERP}(DpF_{
m i+1})$$

where INTERP(x) is the value interpolated from the function $X = D_{avg,j}$ and $Y = dA/dD_j$.

mi micromeritics®

(9) dA/d log(D) pore area [dA/d log(DpF_i), m^2/g]:

$$rac{dA}{d\log(DpF_{
m i})} = {
m INTERP}(DpF_{
m i+1})$$

where *INTERP*(*x*) is the value interpolated from the function $X = D_{avg,j}$ and $Y = dA/d \log(D)_j$.

COMPENDIUM OF VARIABLES

ΔTd	=	thickness of layer desorbed from walls of newly opened pores (Å)
ΔTw	=	thickness of adsorbed layer desorbed during interval (Å)
Α	=	adsorbate property factor; from the BJH Adsorptive Options window
CSA	=	analysis gas molecular cross-sectional area (nm ²), user-entered on the A <i>dsorpt-</i> <i>ive Properties</i> window
CSAa	=	annular cross-sectional area of the desorbed layer (cm ²)
CSAc	=	cross-sectional area of opening of newly opened pores (cm^2)
D _{avg}	=	average pore diameter (Å)
Dp	=	pore (or core) diameter (Å)
F	=	fraction of pores open at both ends; from the <i>BJH Adsorption Report Options</i> win- dow or the <i>BJH Desorption Report Options</i> window
LP	=	length of pore (cm/g)
Ρ	=	relative pressure
Q	=	quantity adsorbed expressed as a volume (cm ³ /g STP)
Rc	=	Kelvin radius (Å) of core
SAw	=	total surface area of walls exposed (cm ² /g)
Tw	=	thickness of remaining adsorbed wall (Å)
Vc	=	volume desorbed from cores of newly opened pores (cm ³ /g)
Vd	=	volume of gas desorbed from walls of previously opened pores (cm 3 /g)
VI	=	liquid equivalent volume of volume adsorbed (cm ³ /g)
V _{mol}	=	liquid molar volume, from the fluid property information

DFT (DENSITY FUNCTIONAL THEORY)

The adsorption isotherm is known to convey a great deal of information about the energetic heterogeneity and geometric topology of the sample under study. The data of physical adsorption have been used for many years as the basis for methods to characterize the surface area and porosity of adsorbents. Real solid surfaces rarely approach ideal uniformity of structure. It is accepted that in general, the surface of even a nonporous material presents areas of greater or lesser attraction for adsorbed molecules.

This energetic heterogeneity greatly affects the shape of the adsorption isotherm with the result that simple theories such as the Langmuir and BET formulas can, at best, give only approximate estimates of surface area. Porous solids virtually are never characterized by a single pore dimension, but instead exhibit a more or less wide distribution of sizes. The observed adsorption isotherm for a typical material is therefore the convolution of an adsorption process with the distribution of one or more properties which affect that process. This was first stated mathematically by Ross and Olivier¹) for the case of surface energy distribution and has become known as the integral equation of adsorption.

THE INTEGRAL EQUATION OF ADSORPTION

In a general form for a single component adsorptive, the integral equation of adsorption can be written as:

$$Q(p) = \int da db dc \dots q(p, a, b, c \dots) \int (a, b, c \dots)$$
(1)

where

Q(p)	 the total quantity adsorbed per unit weight at pressure p,
a,b,c,	 a set of distributed properties,
f(a,b,c,)	 the distribution function of the properties, and
q(p,a,b,c,)	 the kernel function describing the adsorption isotherm on unit surface of material with fixed properties a, b, c,

Equation (1), a Fredholm integral of the first kind, is a member of a class of problems known as illposed, in that there are an infinite number of functional combinations inside the integral that will provide solutions. Even when the kernel function is known, experimental error in the data can make solving for even a single distribution function a difficult task. Solving for multiple distribution functions requires more data than provided by a single adsorption isotherm.

¹) Ross and Olivier, J.P., "On Physical Adsorption," J. Wiley and Sons, New York (1964).

APPLICATION TO SURFACE ENERGY DISTRIBUTION

mi micromeritics®

Under certain conditions, an energetically heterogeneous surface may be characterized by a distribution of adsorptive energies. The conditions are that the sample is not microporous, i.e., that adsorption is taking place on essentially a free surface with no pore filling processes at least to about 0.2 relative pressure. Secondly, that each energetically distinct patch contributes independently to the total adsorption isotherm in proportion to the fraction of the total surface that it represents. This condition is satisfied if the patches are relatively large compared to an adsorptive molecule, or if the energy gradient along the surface is not steep. In mathematical terms, this concept is expressed by the integral equation of adsorption in the following form:

$$Q(p) = \int d\epsilon ~ q(p,\epsilon) \int (\epsilon)$$

(2)

where

- Q(p) = the experimental quantity adsorbed per gram at pressure p,
- $q(p,\epsilon)$ = the quantity adsorbed per unit area at the same pressure, *p*, on an ideal free surface of energy ϵ , and
- $f(\varepsilon)$ = the total area of surface of energy ε in the sample.

The exact form of the energy-dependent term depends on the form of the model isotherms expressed in the kernel function and is provided in the model description.

APPLICATION TO PORE SIZE DISTRIBUTION

Similarly, a sample of porous material may be characterized by its distribution of pore sizes. It is assumed in this case that each pore acts independently. Each pore size present then contributes to the total adsorption isotherm in proportion to the fraction of the total area of the sample that it represents. Mathematically, this relation is expressed by:

$$Q(p) = \int dH \; q(p,H) \int (H)$$

(3)

where

q(p) = the experimental quantity adsorbed at pressure p,
 q(p,H) = the quantity adsorbed per unit area at the same pressure, p, in an ideal pore of size H, and
 f(H) = the total area of pores of size H in the sample.

Numerical values for the kernel functions in the form of model isotherms can be derived from modern statistical mechanics such as density functional theory or molecular simulations, or can be calculated from one of various classical theories based on the Kelvin equation. Several types are found in the models library.

PERFORMING THE DECONVOLUTION

The integrations in equations (2) and (3) are carried out over all surface energies or pore sizes in the model. The functions $q(p,\varepsilon)$ and q(p,H), which we call the kernel functions, are contained in numeric form as model isotherms. Because, in general, there is no analytic solution for equation (1), the problem is best solved in a discrete form; the integral equation for any distributed property Z becomes a summation:

$$Q(p) = \sum_{i} q(p, Z_i) \int (Z_i)$$
(4)

Given a set of model isotherms, q(p,Z), from a model chosen from the models library and an experimental isotherm, Q(p), contained in a sample information file, the software determines the set of positive values f(Z) that most nearly, in a least squares sense, solves equation (4). The distributed property, surface energy or pore size, is then displayed on the *Report Options* window as a selection of tables or graphs.

REGULARIZATION

DFT allows a selectable regularization (also referred to as smoothing) constraint to be applied during the deconvolution process to avoid over-fitting in the case of noisy data or ill-fitting models. The method used is based on co-minimization of the second derivative of the distribution. The relative weight given to this term is determined by the value of the regularization parameter, which is set on the *DFT Pore Size* or *Surface Energy* window and also is shown in the header of reports. The value of the regularization parameter varies from zero (for no second derivative constraint) to ten (indicating a weight equal to minimizing the residuals), or even larger. When the distribution and residuals obtained change little with the value of the regularization parameter, it indicates that the chosen model provides a good representation of the data. Conversely, a large sensitivity to the regularization parameter might indicate inadequate data or a poor choice of model to represent the data.

DOLLIMORE-HEAL PORE VOLUME AND AREA DISTRIBUTION

The calculations for the Dollimore-Heal reports are the same as those for BJH, except for the calculation of average pore diameter and pore length.

PORE DIAMETER

Pore diameter is determined from the Kelvin radius and thickness equation:

$$D_{\mathrm{i}}=2r_{\mathrm{k}}P_{i}+tP_{\mathrm{i}}$$

The average pore diameter is the arithmetic mean of the diameters that bound the interval.

$$\overline{D_{\mathrm{i}}} = \left(rac{D_{i}+D_{\mathrm{i+1}}}{2}
ight)$$

PORE LENGTH

$$egin{aligned} &l_{\mathrm{i}} = rac{A_{\mathrm{p},\mathrm{i}} + 10^8}{\pi \overline{D}_{\mathrm{i}}} \ &A_{\mathrm{p},\mathrm{i}} = rac{4 imes (10^8 \Delta V_{\mathrm{p}})}{\overline{D}_{\mathrm{i}}} \ &\Delta V_{\mathrm{p}} = C
u (D(Q_{\mathrm{i}-1} - Qi) - \Delta t imes 10^8 (A_{\mathrm{p},\mathrm{cum}} - 2\pi \overline{t} l_{\mathrm{i},\mathrm{cum}})) \ &C v = \left(rac{ar{D}_{\mathrm{i}}}{2(ar{\pi}_{\mathrm{k}} + t(P_{\mathrm{i}}) - t(P_{\mathrm{i}+1}))}
ight)^2 \ &ar{t} = rac{ar{D}_{\mathrm{i}}}{2 - ar{ au}_{\mathrm{k}}} \ &ar{ au}_{\mathrm{k}} = rac{(r_{\mathrm{k},\mathrm{i}} + r_{\mathrm{k},\mathrm{i}+1})}{2} \end{aligned}$$

mi micromeritics[®]

where

=	Change in pore volume
=	Pore surface area
=	Summations over the lengths and areas calculated so far
=	Volume correction factor
=	Density conversion factor
=	Average Kelvin radius
=	Average thickness
	= = =

mi micromeritics®

DUBININ-ASTAKHOV

The Dubinin-Astakhov equation is:

$$\log(Q) = \log(Q_0) - \left(rac{RT}{eta E_0}
ight)^N imes \left(\lograc{P_0}{P}
ight)^N$$

where

- β = the affinity coefficient of the analysis gas relative to the P₀ gas, from the *Dubinin* Adsorptive Options window
- E_0 = characteristic energy (*kj/mol*)
- *N* = Astakhov exponent, may be optimized or user entered from the *Dubinin Report Options* window
- *P* = equilibrium pressure
- P_0 = saturation vapor pressure of gas at temperature T
- Q = quantity adsorbed at equilibrium pressure (cm³/g STP)
- Q_0 = the micropore capacity (cm³/g STP)
- R = the gas constant (0.0083144 kj/mol)
- T = analysis bath temperature (K)

For each point designated for Dubinin-Astakhov calculations, the following calculations are done:

$$LV = \log(Q)$$

 $LP = \log(\frac{P_0}{P})^N$

A least-squares fit is performed on the (LP,LV) designated pairs where LP is the independent variable and LV is the dependent variable. If the user selected Yes for the *Optimize Astakhov Exponent* prompt, a systematic search for the optimum value of N is conducted by recalculating the linear regression and selecting the value of N that gives the smallest standard error of the y-intercept. The exponent N is optimized to within 10^{-4.} If the optimum value for N is not found in this range, an exponent of 2 is used. The following are calculated:

- a. Slope (S cm³/g STP)
- b. Y-intercept (YI cm³/g STP)
- c. Error of the slope ($S_{err} \, cm^3/g \, STP$)
- d. Error of the y-intercept ($YI_{err} cm^3/g STP$)
- e. Correlation coefficient
- f. Optimized Astakhov exponent (N)

Using the results of the above calculations, the following can be calculated:

Monolayer Capacity (cm³/g STP):

$$Q_0=10^{YI}$$

Micropore Volume (cm³/g):

$$V_{
m i}=rac{Q_iV_{
m mol}}{22414}$$

where

 V_{mol} = liquid molar volume conversion factor from the fluid property information

Limiting Micropore Volume (cm³/g):

$$V_0=rac{Q_0V_{
m mol}}{22414{
m cm}^3{
m STP}}$$

where

 V_{mol} = liquid molar volume from the fluid property information

Error of Limiting Micropore Volume (cm³/g):

$$V_{0,\mathrm{err}} = W_0 (10 Y I_{\mathrm{err}} - 1.0)$$

Characteristic Energy (KJ/mol):

$$E = rac{2.303(RT)}{eta (2.303 imes S)^{1/N}}$$

mi micromeritics[®]

Modal Equivalent Pore Diameter (Å):

$$D_{
m mode} = 2 \Biggl\{ \left[rac{3N}{3N+1}
ight]^{1/N} imes \left[rac{10^3 nm^3/\dot{A}^3}{eta \cdot E_0}
ight] \Biggr\}^{1/3}$$

where

 β = affinity coefficient of the analysis gas relative to the P₀ gas from the *Dubinin* Adsorptive Options window

Maximum Differential Pore Volume (cm³/g-Å):

This value is also known as frequency of the mode¹).

$$rac{dV}{dD_{
m mode}}{
m Max} = 0.5(3N+1)W_0 igg[rac{3N+1}{3N}igg]^{1/3{
m N}} igg[rac{eta\cdot E_0}{((10^3nm^3)/\dot{A}^3)}igg]^{1/3} {
m exp}igg(-igg[rac{3N+1}{3N}igg]igg)$$

Mean Equivalent Pore Width (Å):

$$D_{ ext{mean}} = 2 imes rac{\left[rac{(10^3 ext{nm}^3) \dot{A}^3}{eta \cdot E_0}
ight]^{1/3}}{ au(rac{3N+1}{3N})}$$

Micropore surface area (m²/g):

$$SDA = 1000 imes 2.0 imes W_0 imes \left[rac{E_0}{k}
ight]^{1/3} imes au igg(rac{3N+1}{3N}igg)$$

 Γ is calculated by a polynomial approximation over the domain $0 \le x \le 1$:

$$egin{aligned} au(x+1) &= 1 + b_1 x + b_2 x^2 + b_3 x^3 + b_4 x^4 + b_5 x^5 + b_6 x^6 + b_7 x^7 + b_8 x^8 \ &+ arepsilon \mathbf{x} |arepsilon \mathbf{x}| \leq 3(10^{-7}) \end{aligned}$$

¹) Ross and Olivier, J.P., "On Physical Adsorption," J. Wiley and Sons, New York (1964).

where

b ₁	=	-0.57719 1652
b ₂	=	0.98820 5891
b ₃	=	-0.89705 6937
b ₄	=	0.91820 6857
b ₅	=	-0.75670 4078
b ₆	=	0.48219 9394
b ₇	=	-0.19352 7818
b ₈	=	0.03586 8343

and where

$$x+1=\left(rac{3N+1}{3N}
ight)$$

Equivalent Pore Diameter (Å):

$$D_{
m i} = 2 igg[rac{- igg(rac{10^3 {
m nm}^3/{
m \AA}^3}{eta E_0} igg)^N}{{
m ln}(W_{
m i}) - {
m ln}(W_0)} igg]^{1/3N}$$

dV/dD Pore Volume (cm³/g-Å):

$$egin{aligned} rac{dV}{dD_{\mathrm{i}}} &= 0.5 imes W_0 imes 3Nigg(rac{10^3\,\mathrm{nm^3/\AA^{3N}}}{eta \cdot E_0}igg)igg(rac{D_{\mathrm{i}}}{2}igg)^{-(3N+1)} \ & imes \expigg[-igg(rac{10^3\,\mathrm{nm^3/\AA^3}}{eta \cdot E_0}igg)^Nigg(rac{D_{\mathrm{i}}}{2}igg)^{-3N}igg] \end{aligned}$$

mi micromeritics®

DUBININ-RADUSHKEVICH

The Dubinin-Radushkevich¹) equation is:

$$\log(Q) = \log(Q_0) - rac{BT^2}{eta} imes \left[\log rac{P_0}{P}
ight]^2$$

where

- β = the affinity coefficient of analysis gas relative to P₀ gas (for this application β is taken to be 1)
- B = a constant
- P_0 = saturation vapor pressure of gas at temperature T
- P = equilibrium pressure
- Q = quantity adsorbed at equilibrium pressure (cm³/g STP)
- Q_0 = the micropore capacity (cm³/g STP)
- T = analysis bath temperature (K), from the P_0 and Temperature Options window

For each point designated for Dubinin-Radushkevich calculations, the following calculations are done:

$$LV = \log(Q)$$

$$LP = \log(\frac{P_0}{P})^2$$

The intercept, log(Vo) can be found by performing a least-squares fit on the (LP,LV) designated pairs where LP is the independent variable and LV is the dependent variable. Assuming the adsorption of gas is restricted to a monolayer, Vo is the monolayer capacity. Based on this assumption, the following are calculated:

¹⁾ Dubinin, M., *Carbon 21*, 359 (1983); Dubinin, M., *Progress in Surface and Membrane Science 9*, 1, Academic Press, New York (1975); Dubinin, M. and Astakhov, V., *Adv. Chem. Ser. 102*, 69 (1971); Lamond, T. and Marsh, H., *Carbon 1*, 281, 293 (1964); Medek, J., *Fuel 56*, 131 (1977); Polanyi, M., *Trans. Faraday Soc. 28*, 316 (1932); Radushkevich, L., *Zh. fiz. Kemi. 33*, 2202 (1949); Stoeckli, H., *et al, Carbon 27*, 125 (1989).

- a. Slope (S cm³/g STP)
- b. Y-intercept (Y/ cm^{3/}g STP)
- c. Error of the slope (S_{err} cm³/g STP)
- d. Error of the y-intercept (YI_{err} cm³/g STP)
- e. Correlation coefficient

Using the results of the above calculations, the following can be calculated:

Monolayer Capacity (cm³/g STP):

$$Q_0=10^{YI}$$

Error of Monolayer Capacity (cm³/g STP):

$$Q_{0,\mathrm{err}} = Q_0 (10^{YI,\mathrm{err}} - 1.0)$$

Micropore surface area (m²/g):

$$SDP = rac{\sigma Q_0 N_A}{22414 {
m cm}^3 \left(rac{10^{18} {
m nm}^2}{{
m m}^2}
ight)}$$

where

 σ = molecular cross sectional area of gas (nm²) from the *Adsorptive Properties* window

EQUILIBRATION

Equilibration is reached when the pressure change per equilibration time interval (first derivative) is less than 0.01% of the average pressure during the interval. Both the first derivative and average pressure are calculated using the Savitzky-Golay¹) convolution method for polynomial functions. The following equations are those used to compute weighted average and first derivative, respectively, for the 6th point of an 11-point window:

$$P_{\text{avg}} = \frac{-36(P_{11}+P_1)+9(P_{10}+P_2)+44(P_9+P_3)+69(P_3+P_4)+84(P_7+P_5)+89(P_6)}{429}$$

 $P_{\rm chg} = \frac{5(P_{11}-P_1) + 4(P_{10}-P_2) + 3(P_9-P_3) + 2(P_8-P_4) + (P_7-P_5)}{110}$

$$P_{ ext{pcp},i} = 100\% rac{P_{ ext{chg}}}{P_{ ext{avg}}}$$

pressure change per equilibration time interval

where the numerical constants are from the Savitzky-Golay convolution arrays, and

P _{avg}	= average pressure
P _{chg}	= change in pressure
P _{pcp,i}	= percent change per interval
Pi	= i th pressure reading taken at equilibrium intervals

If a non-zero value that is too small is entered for the maximum equilibration time, the points are collected before equilibration is reached.

If P_{avg} is greater than 0.995 times the current P_0 , equilibration will not take place until the *Minimum equilibration delay for P/P*₀ 0.995 has expired, in addition to the standard equilibration criteria.

¹) Savitzky, A. and Golay, M.J.E., Anal. Chem. 36, 1627 (1964).

F-RATIO METHOD

The *f*-Ratio is the quantity adsorbed divided by the quantity adsorbed in a reference isotherm at the same pressure.

$$f_i = rac{Q_i}{Q_{ ext{ref}}P_i}$$

The reference quantity adsorbed is found by spline interpolation of the reference isotherm.

FREUNDLICH ISOTHERM

The Freundlich isotherm has the form:

$$rac{Q}{Q_{ extsf{m}}} = CP^{rac{1}{ extsf{m}}}$$

where

C = temperature-dependent constant

m = temperature-dependent constant

P = equilibrated collected pressure measured by gauge at temp T_{amb}

Q = quantity of gas adsorbed

 Q_m = quantity of gas in a monolayer

The pressure is absolute; typically, m > 1. In terms of quantity adsorbed,

 $Q = Q_{\rm m} C P^{rac{1}{{
m m}}}$

Taking the log of both sides yields

 $\log Q = \log Q_{\rm m}C + \frac{1}{{
m m}}\log P$

HEAT OF ADSORPTION FOR PHYSICAL ADSORPTION

The adsorption isostere is represented by

$$\ln\left(rac{P}{P_0}
ight) = rac{q_i}{RT} + C$$

where

- q_i = isosteric heat of adsorption
- C = unknown constant

R = gas constant

T = analysis temperature

The isosteric heat of adsorption, q_i , is determined by finding the slope of *In* (*P*/*P*0) as a function of 1/RT for a set of isotherms measured at different temperatures.

Horvath-Kawazoe

A relative pressure lower limit is determined such that L- d_0 never equals zero. All pressure points less than this limit are discarded. For each collected relative pressure point, values of L are chosen in an iterative manner, and the relative pressure (P/P_0) determined by solving one of the following equations:

- Slit Pore Geometry (original Horvath-Kawazoe)
- Cylinder Pore Geometry (Saito / Foley)
- Sphere Pore Geometry (Cheng / Yang)

SLIT PORE GEOMETRY (ORIGINAL HORVATH-KAWAZOE)

When using the original Horvath-Kawazoe¹) method, the following equation is solved for each value of *P*. The value of *L* is determined when the solved-for relative pressure is within 0.1% of the collected absolute pressure:

$$\ln \frac{P}{P_0} = \frac{N_A}{RT} \frac{IP \times 10^{32}}{\sigma^4 L - 2d_0} \frac{\sigma_4}{3L - {d_0}^3} - \frac{\sigma^{10}}{9L - {d_0}^9} - \frac{\sigma^4}{3{d_0}^3} + \frac{\sigma^{10}}{9{d_0}^9}$$

where

$$10^{32}$$
 = the number of cm⁴ that are equal to Å⁴

$$Z_{\rm s}+Z_{\rm A}$$

$$\sigma$$
 = gas solid nuclear separation at zero interaction energy (Å), $\frac{1}{2}$

$$d_0 = \frac{D_A + D_s}{2}$$

where

- D_A = molecular diameter (Å) from the *Horvath-Kawazoe Physical Properties* window
- D_S = diameter of sample atom (Å) from the *Horvath-Kawazoe Physical Properties* window

IP = interaction parameter (erg-cm⁴) from the *Horvath-Kawazoe Report Options* window

 N_A = Avogadro Constant (6.02214129 × 10²³ mol⁻¹)

L = pore width (nucleus to nucleus) (Å)

¹) Horvath, G. and Kawazoe, K., *J. Chem. Eng. Japan 16(6)*, 470 (1983).

- P = equilibrium pressure
- P₀ = saturation pressure
- R = gas constant (8.31441 × 10^7 erg/mol K)
- T = analysis bath temperature (K), from an entered or calculated value on the P_0 and *Temperature Options* window

where:

- Z_s = sample equilibrium diameter at zero interaction energy (Å) from the *Horvath-Kawazoe Physical Properties* window
- Z_A = zero interaction energy diameter from the *Horvath-Kawazoe Physical Properties* window

CYLINDER PORE GEOMETRY (SAITO/FOLEY)

When using the Saito/Foley¹) method, the following equation is solved for each value of *P*. The value of *L* is determined when the solved-for relative pressure is within 0.1% of the collected absolute pressure:

$$egin{aligned} &\ln\left(rac{P}{P_0}
ight) = rac{3}{4}rac{\pi N_A}{RT} imes rac{IP imes 10^{32}}{{d_0}^4} imes \sum_{k=0}^\infty \left[rac{1}{ ext{k}+1} igg(1-rac{d_0}{r_ ext{p}}igg)^{2k} imes \ &\left\{rac{21}{32}lpha_ ext{k} igg(rac{d_0}{r_ ext{p}}igg)^{10} - eta_ ext{k} igg(rac{d_0}{r_ ext{p}}igg)^4
ight\}
ight] \end{aligned}$$

where

 10^{32} = the number of Å⁴ that are equal to 1 cm⁴

$$\begin{array}{lll} \beta_{k} & = & \left(\frac{-1.5-k}{k}\right)^{2}\beta_{k-1}, \beta_{0} = 1.0\\ \\ \alpha_{k} & = & \left(\frac{-4.5-k}{k}\right)^{2}\alpha_{k-1}, \alpha_{0} = 1.0 \end{array}$$

¹) Saito, A. and Foley, H. C., *AIChE Journal 37(3)*, 429 (1991).

$$d_0 = \frac{D_A D_s}{2}$$

where

- D_A = molecular diameter (Å) from the *Horvath-Kawazoe Physical Properties* window
- D_S = diameter of sample atom (Å) from the *Horvath-Kawazoe Physical Properties* window
- $IP = \frac{\text{interaction parameter (10^{-43} erg-cm^4) from the Horvath-Kawazoe Report Options window}}{\text{dow}}$
- N_A = Avogadro Constant (6.02214129 × 10²³ mol⁻¹)
- L = pore width (nucleus to nucleus) (Å)
- P = equilibrium pressure
- P_o = saturation pressure
- R = gas constant (8.31441 × 10^7 erg/mol K)
- r_p = radius of the cylindrical pore, $\frac{L}{2}$
- T = analysis bath temperature (K), from an entered or calculated value on the P_0 and Temperature Options window

SPHERE PORE GEOMETRY (CHENG/YANG)

When using the Cheng / Yang¹) method, the following equation is solved for each value of P. The value of L is determined when the solved-for relative pressure is within 0.1% of the collected absolute pressure:

$$\ln\left(\frac{P}{P_0}\right) = \frac{6N_1\epsilon^*_{12} + N_2\epsilon^*_{22}L^3 \times 10^{32}}{RTL - d_0^3} \left[-\left(\frac{d_0}{L}\right)^6 \left(\frac{T_1}{12} + \frac{T_2}{8}\right) + \left(\frac{d_0}{L}\right)^{12} \left(\frac{T_3}{90} + \frac{T_4}{80}\right) \right]$$

where

 10^{32} = the number of cm⁴ that are equal to Å⁴

$$\epsilon^{*}_{12} = \frac{\dot{A}_{S}}{4d_{s}^{6}}, \text{where}\dot{A}_{s} = \frac{6(mc^{2})\alpha_{S}\alpha_{A}}{\frac{\alpha_{S}}{X_{S}} + \frac{\alpha_{A}}{X_{A}}}$$

 $\epsilon^{*}_{22} = \frac{\dot{A}_{A}}{4D_{A}^{6}}, \text{where}\dot{A}_{A} = \frac{3(mc^{2})(\alpha_{A})(\chi_{A})}{2}$

$$d_0 = \frac{D_A + D_S}{2}$$

where

D_A = molecular diameter (Å) from the *Horvath-Kawazoe Physical Properties* window

D_S = diameter of sample atom (Å) from the *Horvath-Kawazoe Physical Properties* window

- L = pore width (nucleus to nucleus) (Å)
- $N_1 = 4\pi L_2 N_S$, where N_s = number of sample atoms/cm² at monolayer

N₂ = $4\pi (L-d_0)^2 N_A$, where N_S = number of gas molecules/cm²

- P = equilibrium pressure
- P_0 = saturation pressure
- R = gas constant (8.31441 × 10^7 erg/mol K)
- T = $analysis bath temperature (K), from an entered or calculated value on the <math>P_0$ and Temperature Options window

¹) Cheng, Linda S. and Yang, Ralph T., *Chemical Engineering Science* 49(16), 2599-2609 (1994).

$$T_1 = \frac{1}{(1-S)^3} - \frac{1}{(1+S)^3}$$
$$T_2 = \frac{1}{(1+S)^2} - \frac{1}{(1-S)^2}$$
$$T_3 = \frac{1}{(1-S)^9} - \frac{1}{(1+S)^9}$$
$$T_4 = \frac{1}{(1+S)^8} - \frac{1}{(1-S)^8}$$
where $S = \frac{L-d_0}{L}$

CHENG/YANG CORRECTION

This factor corrects for the nonlinearity of the isotherm. It adds an additional term to the equations for the different geometrics:

$$\ln\left(rac{P}{P_0}
ight) = G(L) - \left[1 - rac{1}{ heta} \ln\left(rac{1}{1- heta}
ight)
ight]$$

where

G(L) = one of the Horvath-Kawazoe equations given above

 θ = degree of void filling; θ is estimated by first computing the monolayer capacity (Q_m) with the Langmuir equation over the range of data points from relative pressure 0.02 to 0.2 or the maximum relative pressure included in the Horvath-Kawazoe analysis. θ is computed as the quantity adsorbed over Q_m .

INTERACTION PARAMETER

The interaction parameter (IP) results from the following calculations:

The Kirkwood-Muller dispersion coefficients

$$egin{aligned} A_{\mathrm{S}} &= rac{6mc^2lpha_{\mathrm{S}}lpha_{\mathrm{A}}}{rac{lpha_{\mathrm{S}}}{\chi_{\mathrm{S}}}+rac{lpha_{\mathrm{A}}}{\chi_{\mathrm{A}}}} \ A_{\mathrm{A}} &= rac{3mc^2lpha_{\mathrm{A}}\chi_{\mathrm{A}}}{2} \end{aligned}$$

where

- α_A = polarizability of gas molecule (cm³)
- $\alpha_{\rm S}$ = polarizability of sample atoms (cm³)
- mc^2 = kinetic energy of electron (0.8183 × 10⁻⁶ erg)
- X_A = diamagnetic susceptibility of gas molecule (cm³)

$$IP = (N_{\mathrm{M}}A_{\mathrm{A}}) + (N_{\mathrm{S}}A_{\mathrm{S}})$$

where:

- N_{M} = number of gas molecules/cm² at monolayer from the *Horvath-Kawazoe Physical Properties* window
- N_S = number of sample atoms/cm² from the *Horvath-Kawazoe Physical Properties* window
- X_s = diamagnetic susceptibility of sample atom (cm³)

See Interaction Parameter Components Table on page 34 for recommended values.

ADDITIONAL CALCULATIONS

Based on the previous calculations, the following can be calculated:

Adjusted Pore Width (Å):

(Shell to Shell)

 $AL_{
m i}=L_{
m i}-D_{
m S}$

Cumulative Pore Volume (cm³/g):

$$V_{ ext{cum,i}} = rac{Q_{ ext{i}} V_{ ext{mol}}}{22414 ext{cm}^3 ext{STP}}$$

where

 V_{mol} = liquid molar volume from the fluid property information

dV/dD Pore Volume (cm³/g-Å):

$$rac{dV}{dD_{
m i}} = rac{V_{
m cum,i}-V_{
m cum,i-1}}{AL_{
m i}-AL_{
m i-1}}$$

Median Pore Width (Å):

$$egin{aligned} V_{ ext{half}} &= rac{V_{ ext{cum,n}}}{2} \ D_{ ext{med}} &= \expigg(\ln D_{ ext{l}} + (\ln D_{ ext{g}} - \ln D_{ ext{l}}) rac{\ln V_{ ext{half}} - \ln V_l}{\ln V_g - \ln V_l}igg) \end{aligned}$$

where

DI	=	pore width (L_i) that corresponds to $V_{\rm I}$
Dg	=	pore width (L_i) that corresponds to V_{g}
V _{cum,n}	=	total cumulative pore volume ($V_{\text{cum},i}$) for points designated for Horvath-Kawazoe calculations
V_g	=	cumulative pore volume ($V_{cum,i}$) for first point greater than V_{half}
V _{half}	=	50% of total cumulative pore volume
VI	=	cumulative pore volume ($V_{cum,i}$) for first point less than V_{half}

Gas	Bath Temperature (K)	Sample Type	Interaction Parameter Calculated Value *
Argon	87.3	Carbon (Ross/Olivier value) Carbon (Horvath/Kawazoe value) Zeolite	2.61 5.89 3.19
Carbon Dioxide	298.15	Carbon (Ross/Olivier value) Carbon (Horvath/Kawazoe value) Zeolite	4.20 9.20 5.08
	273.15	Carbon (Ross/Olivier value) Carbon (Horvath/Kawazoe value) Zeolite	4.34 9.35 5.22
	194.65	Carbon (Ross/Olivier value) Carbon (Horvath/Kawazoe value) Zeolite	4.72 9.72 5.60
Nitrogen	77.15	Carbon (Ross/Olivier value) Carbon (Horvath/Kawazoe value) Zeolite	2.84 6.53 3.49
following fiel	d:	entered in the <i>Horvath-Kawazoe Report Op</i> Ited value) × 10 ⁻⁴³ erg-cm ⁴	<i>tions</i> window in the

The following values were used to calculate the values in the previous table.

Carbon-Graphite

Zeolite

DS	=	3.40	DS	=	3.04
NS	=	3.845 × 10 ¹⁵	NS	=	3.75 × 10 ¹⁵
X _S	=	1.05 × 10 ⁻²⁹ (Ross/Olivier)	X _S	=	1.94 × 10 ⁻²⁹
		13.5 × 10 ⁻²⁹ (Horvath/Kawazoe, implicit)	a_{s}	=	0.85 × 10 ⁻²⁴
α _s	=	1.02 × 10 ⁻²⁴			

Nitrogen

Argon

D_A	=	3.00	D _A	=	2.95
N _A	=	6.71 × 10 ¹⁴	N _A	=	7.608 × 10 ¹⁴
X _A	=	3.6 × 10 ⁻²⁹	X _A	=	3.22 × 10 ⁻²⁹
α _A	=	1.76 ×10 ⁻²⁴	α _A	=	1.63 ×10 ⁻²⁴

Carbon Dioxide

D_A	=	3.23
N _A	=	4.567 × 10 ¹⁴ (25 ℃)
		5.45 × 10 ¹⁴ (0 ℃)
		7.697 × 10 ¹⁴ (-78 ℃)
X _A	=	5.0 × 10 ⁻²⁹
α _A	=	2.7 × 10 ⁻²⁴

 D_A values are from van der Waal's constant.

 N_A values are from liquid densities.

X and α values are derived from data found in Ross and Olivier¹).

¹) Ross and Olivier, J.P., "On Physical Adsorption," J. Wiley and Sons, New York (1964)

The physical parameters referenced in Saito/Foley are:

Alumino	phos	phate		Aluminosilicate		
C) _S	=	2.60	D _S	=	2.76
Ν	۱ _S	=	1.48 × 10 ¹⁵	N _S	=	1.31 × 10 ¹⁵
X	(_S	=	1.3 × 10 ⁻²⁹	X _S	=	1.3 × 10 ⁻²⁹
C	۱ _s	=	2.5 × 10 ⁻²⁴	α _s	=	2.5 × 10 ⁻²⁴

LANGMUIR SURFACE AREA

For each point designated for surface area calculations, the Langmuir¹) transformation is calculated as:

$$L=rac{P_{
m rel}}{N_{
m ads}}$$

where *L* is in units of g/cm^3 STP.

A least-squares fit is performed on the (P_{rel} , L) designated pairs where P_{rel} is the independent variable and L is the dependent variable. The following are calculated:

- a. Slope (S g/cm³ STP)
- b. Y-intercept (Y_{int} g/cm³ STP)
- c. Error of the slope (S_{err} g/cm³ STP)
- d. Error of the y-intercept (YI_{err} g/cm³ STP)
- e. Correlation coefficient

Using the results of the above calculations, the following can be calculated:

Langmuir Surface Area (m²/g):

$$SA_{
m LAN} = rac{{
m CSA_{ imes}}N_A}{(22414 {
m cm}^3 {
m STP})(rac{10^{18}nm^2}{m^2})S}$$

where

- *CSA* = analysis gas molecular cross-sectional area (nm²), user-entered on the *Adsorptive Properties* window
- N_A = Avogadro constant (6.02214129 × 10²³ mol⁻¹)

Quantity of the Monolayer (cm³/g STP):

$$Q_m = rac{1}{s}$$

¹) Langmuir, I., J. *Am. Chem. Soc. 38*, 2267 (1916); J. Am. Chem. Soc. 40, 1361 (1918); Phys. Rev 8, 149 (1916).



Langmuir b Value:

$$b = Y_{\mathrm{int}} V_m$$

Error of the Langmuir Surface Area (m²/g):

$$LAN_{
m err} = rac{SA_{
m Lan}S_{
m err}}{S}$$

MP-METHOD

With the (t_i, Q_i) data pairs¹), the Akima semi-spline interpolation method is used to interpolate quantity adsorbed values based on thickness values that are evenly spaced 0.2 angstrom apart starting at the first outlier point. Outliers are defined as those points that have the maximum instantaneous slope within an iteratively shrinking subset of all points. The remaining pore surface area calculation result is the slope of the line defined by two consecutive interpolated points. The slopes of each pair of consecutive points from the origin to the last point must be monotonically decreasing and non-negative. With the interpolated points set the following can be calculated:

Average pore hydraulic radius (Å):

$$R_{\mathrm{i}} = rac{t_{\mathrm{i}} + t_{\mathrm{i-1}}}{2}$$

Remaining pore surface area for the ith point (m²/g):

$$S_{
m i} = rac{Q_{
m i} - Q_{
m i-1}}{t_{
m i} - t_{
m i-1}} rac{V_{
m mol}}{22414 {
m cm}^3 {
m STP}} imes 10^4$$

where

 10^4 = unit conversions V_{mol} = liquid molar volume from the fluid property information

Incremental pore surface area occluded for the ith point (m²/g):

$$S_{\mathrm{inc},i} = S_{\mathrm{i}-1} - S_{\mathrm{i}}$$

Cumulative pore surface area occluded for the ith point (m²/g):

 $S_{ ext{cum}_{ ext{i}}} = S_{ ext{inc}, ext{i}} + S_{ ext{inc}, ext{i}-1} + \ldots + S_{ ext{inc}, ext{i}}$

dA/dR pore surface area for the i^{th} point (m²/g-Å):

$$rac{dA}{dR_{
m i}}=rac{S_{
m inc_{
m i}}}{t_{
m i}-t_{
m i-1}}$$

¹) Mikhail, R., Brunauer, S. and Bodor, E., J. Colloid and Interface Sci. 24, 45 (1968).

Incremental pore volume occluded for the ith point (cm³/g):

$$V_{
m inc,i}=S_{
m inc,i}R_i imes 10^{-4}$$

Cumulative pore volume occluded for the ith point (cm³/g):

$$V_{ ext{cum}, ext{i}} = V_{ ext{inc}, ext{i}-1} + V + \ldots + V_{ ext{inc}, ext{i}}$$

dV/dR pore volume for the i^{th} point (cm³/g-Å):

$$rac{dV}{dR_{
m i}}=rac{V_{{
m inc},i}}{t_{
m i}-t_{
m i-1}}$$

QUANTITY ADSORBED

FREE SPACE - MEASURED

Measured free-space volumes are calculated using the following equations:

$$egin{aligned} V_{fw} &= \left(rac{P_1}{P_2}V_{man1} - V_{man2}
ight)rac{T_{STD}}{T_{man}} \ V_{fc} &= V_{fw}\left(rac{P_2}{P_3}
ight) \end{aligned}$$

$$V_{ ext{bath}} = rac{V_{fc} - V_{fw}}{1 - rac{T_{ ext{bath}}}{T_{ ext{amb}}}}$$

FREE SPACE - CALCULATED

The calculated free space is determined by subtracting the gas capacity of the volume occupied by the sample from the measured free space of the empty tube.

$$egin{aligned} V_{fw} &= V_{wb} + V_p - V_s \left(rac{T_{ ext{STD}}}{T_{ ext{amb}}}
ight) \ V_{fc} &= V_{cb} + V_p - V_s \left(rac{T_{ ext{STD}}}{T_{ ext{bath}}}
ight) \end{aligned}$$

COMPENDIUM OF VARIABLES

P ₁	=	system manifold pressure before dosing onto sample
P ₂	=	system manifold pressure after dosing onto sample
P ₃	=	sample pressure after raising Dewar and equilibrating with helium
T _{amb}	=	approximate room temperature (298 K)
T _{bath}	=	analysis bath temperature (K)
T _{man}	=	system manifold temperature (K)
T _{STD}	=	standard temperature (273.15 K)
V _{bath}	=	volume of cold free space at analysis bath temperature (cm ³ at standard temperature)
V _{cb}	=	volume of cold free space of the empty tube (cm^3 at standard temperature)
V _{fc}	=	volume of cold free space (cm ³ at standard temperature)
V _{fw}	=	volume of warm free space (cm ³ at standard temperature)
V _{man1}	=	manifold volume (cm ³) before dosing onto sample. May differ from V_{man2} because it includes deflection if 10 mmHg transducer is present.
V _{man2}	=	manifold volume (cm 3) after dosing onto sample. May differ from V _{man1}
Vp	=	port volume
Vs	=	sample volume
V_{wb}	=	volume of warm free space of the empty tube (cm 3 at standard temperature)

RELATIVE PRESSURE

If P_0 is measured in the P_0 tube, the current pressure is measured in the P_0 tube when each point is taken, and used to calculate relative pressure for that point:

$$P_{\mathrm{rel}} = rac{P}{P_o}$$

SATURATION PRESSURE

Saturation pressure (P_0) is selected on the P_0 and Temperature Options window. It may be entered or measured in the P_0 tube.

- If entered, P₀ equals user-entered value.
- If measured, P₀ equals equilibrated pressure reading after saturatingP₀ tube or the sample tube.
 P₀ is measured for each data point.

SPC REPORT VARIABLES

REGRESSION CHART VARIABLES

The line of best fit for the Regression Chart is calculated by the usual least squares method. ¹) If there is only a single point or all N points have the same x-value, there can be no line of best fit in the standard form.

$$\bar{\chi} = \frac{\sum x_i}{N}$$
$$\bar{y} = \frac{\sum y_i}{N}$$
Slope = $\frac{\sum (x_i - \bar{\chi})(y_i - \bar{y})}{\sum (x_i - \bar{\chi})^2}$

$$ext{Intercept} = ar{y} - ext{Slope} \cdot ar{\chi}$$

The coefficient of correlation for this line is also calculated in the usual way. 2)

$$\sigma_{\chi} = \sqrt{rac{\sum (\chi_4 - ar{\chi})^2}{N}}$$
 $\sigma_y = \sqrt{rac{\sum (y_i - ar{y})^2}{N}}$

$$\mathrm{Cov}(x,y) = rac{\sum (\chi_i - ar{\chi})(y_i - ar{y})}{N}$$

 $ext{CorrelationCoeff} = rac{ ext{Cov}(x,y)}{\sigma_\chi \sigma_y}$

¹) *BASIC Scientific Subroutines Vol II*, by F.R. Ruckdeschel, Copyright 1981 BYTE Publications/McGraw Hill, p. 16.

²) *Mathematical Handbook for Scientists and Engineers*, G.A. Korn and T.M. Korn, McGraw Hill, Sec. 18.4. (1968)

CONTROL CHART VARIABLES

 $Mean = \frac{\sum y_i}{N}$

 $ext{StandardDeviation} = \sqrt{rac{\sum (y_i - ext{Mean})^2}{N-1}}$

 $C.V. = \frac{StdDev}{Mean}$

 $+n\sigma = \text{Mean} + n \cdot \text{StandardDeviation}$

 $-n\sigma = \text{Mean} - n \cdot \text{StandardDeviation}$

SUMMARY REPORT

The following calculations and the results of previous calculations (as noted) are used to generate the summary report:

a. Single-point Surface Area (m²/g)

$$S_{1PT} = rac{[Q(1-P)] imes CSA(6.023 imes 10^{23})}{22414 ext{cm}^2 imes ext{STP}igg(rac{10^{18} ext{nm}^2}{ ext{m}^2}igg)}$$

where

P = pressure closest to 0.3 of the relative pressure points designated for surface area calculations

Q = quantity adsorbed corresponding to P

- b. Multi-point Surface Area.See **BET Surface Area on page 2**.
- c. Langmuir Surface Area. See Langmuir Surface Area on page 37
- d. t-Plot Micropore Surface Area. See *t-Plot on page 49*.
- e. t-Plot External Surface Area. See t-Plot on page 49.
- f. BJH Cumulative Adsorption
- g. BJH Cumulative Desorption
- h. Adsorption Total Pore Volume
- i. Desorption Total Pore Volume
- j. t-Plot Micropore Pore Volume. See <u>t-Plot on page 49</u>.
- k. Freundlich. See Freundlich Isotherm on page 25.
- I. Temkin. See <u>Temkin Isotherm on page 50</u>.
- m. Alpha-S. See .
- n. DFT Pore Size and DFT Surface Energy. See <u>DFT (Density Functional Theory) on</u> page 12.
- o. Nanoparticle Size

$$d=rac{6 imes 10^4}{A
ho}$$

where

- ρ = sample density
- A = BET surface area
- d = side length (for cubic particles or diameter (for spherical particles)
- p. Dubinin-Astakhov Micropore Surface Area. See Dubinin-Astakhov on page 17.
- q. Dubinin-Astakhov Micropore Volume. See <u>Dubinin-Astakhov on page 17</u>.

- r. Dubinin-Radushkevich Micropore Surface Area. See <u>Dubinin-Radushkevich on page 21</u>.
- s. Dubinin-Radushkevich Monolayer Capacity. See <u>Dubinin-Radushkevich on page 21</u>.
- t. MP-Method Cumulative Surface Area of Pores

 $S_{
m total} = S_{
m cum,i}$

See <u>MP-Method on page 39</u> for the last collected data point used in the MP-method Calculations, and the range of hydraulic pore radii over which the cumulative surface area was computed.

u. MP-Method Cumulative Pore Volume of Pores $V_{
m total} = V_{
m cum.i}$

 $v_{\rm total} = v_{\rm cum,i}$

See <u>MP-Method on page 39</u> for the last collected data point used in the MP-method calculations, and the range of hydraulic pore radii over which the cumulative pore volume was computed.

v. Average Pore Hydraulic Radius (Å)

$$ar{r} = rac{V_{ ext{total}}}{S_{ ext{total}}} imes 10^4$$

w. Horvath-Kawazoe. See Horvath-Kawazoe on page 27.

т**-Р**LОТ

A least-squares analysis fit is performed on the $(t_i, N_{ads,i})$ data pairs where t_i is the independent variable and $N_{ads,i}$ is the dependent variable. Only the values of t_i between t_{min} and t_{max} , the minimum and maximum thickness, are used. The following are calculated:

- a. Slope (S cm³/g-Å STP)
- b. Y-intercept (Y_{int} cm³/g STP)
- c. Error of the slope (S_{err} cm³/g-Å STP)
- d. Error of the Y-intercept (YI_{err} cm³/g STP)
- e. Correlation coefficient

Using the results of the above calculations, the following can be calculated:

External Surface Area (m²/g):

 $rac{SV_{
m mol}}{F imes 22414 {
m cm}^3 STP} imes 10^4$

where

10⁴ = unit conversions

- *F* = surface area correction factor, user-entered on the *t-Plot Report Options* window
- V_{mol} = liquid molar volume, from the fluid property information

Micropore Surface Area (m²/g):

 $SA_{\mu
ho}=SA_{
m total}+SA_{
m ext}$

where SA_{total} is the BET surface area if the user enabled the BET report exclusively, or Langmuir surface area if the user enabled the Langmuir report exclusively. If neither report has been selected, SA_{total} is the BET surface area value calculated using a set of default parameters.

Micropore Volume (cm³ liquid/g):

 $\frac{Y_{\rm int}V_{\rm mol}}{22414{\rm cm}^3{\rm STP}}$

Temkin Isotherm

The Temkin isotherm has the form:

$$rac{Q}{Q_{ extrm{m}}} = rac{RT}{q_{0}lpha} extrm{ln}(A_{0}P)$$

where

 A_0 = adjustable constant

- *a* = adjustable constant
- P = equilibrium pressure measured by gauge at temp T_{amb}
- q_0 = the differential heat of adsorption at zero surface coverage

Q = quantity of gas adsorbed

- Q_m = quantity of gas in a monolayer
- R = molar gas constant 8.31441 $\times 10^{-3} \frac{kJ}{mol K}$

T = bath temperature

In terms of quantity adsorbed

$$Q=rac{RTQ_m}{q_0lpha}[\ln A_0+\ln P]$$

Thus, the plot of the natural log of absolute pressure vs. quantity adsorbed yields a straight line with slope $\frac{RTQ_m}{q_0}$ and intercept $\ln A_0 \frac{RTQ_m}{q_0 \alpha}$

mi micromeritics[®]

THERMAL TRANSPIRATION CORRECTION

During data reduction, thermal transpiration correction is applied to the data if the user selected *Apply thermal transpiration correction* from the *Report Options* window. Starting with the first collected pressure, the following calculations are performed until the pressure ratio (PC/P) is greater than or equal to 0.99.

$$egin{aligned} Y &= \left(rac{P imes ext{SD} imes ext{MD}^2}{2.33 imes T}
ight) 10^3 \ \mu &= rac{(1+G)Y}{(1+H)Y} \ F &= rac{1}{lpha Y^2 + eta Y + \mu} \ P &= \left(1 - F\left(1 - \sqrt{rac{T_{ ext{bath}}}{T_{ ext{amb}}}}
ight)
ight) \end{aligned}$$

where

α	=	Weber's coefficient, 0.033
β	=	Weber's coefficient, 0.245
<i>F</i> ,Υ, μ	=	intermediate values for subsequent calculations
G	=	Weber's coefficient, 2.5
Н	=	Weber's coefficient, 2
MD	=	thermal transpiration hard sphere diameter of gas (Å), from the Adsorptive Prop- erties window
Ρ	=	equilibrated collected pressure measured by gauge at temp $T_{ m amb}$
SD	=	inside diameter of sample tube (mm) from the <i>Report Options</i> window
Т	=	average temperature $rac{T_{ ext{bath}}+T}{2}$
T _{amb}	=	room temperature (298 K)
T _{bath}	=	analysis bath temperature (K), from the P_0 and Temperature Options window

THICKNESS CURVE

For each point designated, the following parameters are used in thickness curve calculations:

parameter #1

- C₂ = parameter #2
- C₃ = parameter #3
- $P_{rel, i}$ = relative pressure for the ith point (mmHg)
- t_i = thickness for ith point

REFERENCE

Interpolated from table.

KRUK-JARONIEC-SAYARI

$$t_{\mathrm{i}} = \left(rac{C_2}{C_2 - \mathrm{log}(P_{\mathrm{rel},i})}
ight)^{c^3}$$

HALSEY

$$t_{
m i} = C_1 igg[rac{C_2}{\ln(P_{
m rel,i})} igg]^{C_3}$$
Halsey¹)

HARKINS AND JURA

$$t_{\mathrm{i}} = \left[rac{C_1}{C_2 - \mathrm{log}(P_{\mathrm{rel},\mathrm{i}})}
ight]^{C_3}$$

Harkins and Jura²)

¹) Halsey, G., J. *Chem. Phys. 16*, 931-937 (1948).

²) Harkins, W.D. and Jura, G., J. *Chem. Phys.* 11, 431 (1943).

BROEKOFF-DE BOER

 $\log(P_{\mathrm{rel},\mathrm{i}}) = rac{C_1}{t^2,\mathrm{i}} + C_2 \exp(c_3 t_\mathrm{i})$

CARBON BLACK STSA

 $t_\mathrm{i} = C_1(P_\mathrm{rel,i})^2 + C_2(P_\mathrm{rel,i}) + C_3$